

The Delta® PetroCut™ Geological Cutter is a manual, 4Hp (3.0kW) cutter using a 10" (254mm) Continuous Rim Diamond Blade for sectioning of rock, minerals, ceramics, coal, glass, concrete, refractory and many geological samples up to 3.75" (95mm) in size. Features include a Rock Clamp Assembly with an X-Y table bed and cross-feed table for serial cutting. Clamps can be positioned at various points on the table bed. Other specimen vises are available for holding irregular shapes.

The Delta® PetroCut™ Geological Cutter includes a clean-out hose for quick debris removal. Two jets optimize presentation of the coolant to the wheel and specimen. Safety is maintained with an emergency stop, hood lock and a mechanical brake which is engaged when the power is turned off. The heavy duty aluminum cast base provides stability and corrosion resistance. Operation instructions are included.

10-2156-260 Delta® PetroCut™, 230 VAC, 60Hz, 3 Phase

10-2156-460 Delta® PetroCut™, 460 VAC, 60Hz, 3 Phase

10-2156-250 Delta® PetroCut™, 220 VAC, 50Hz, 3 Phase

10-2156-400 Delta® PetroCut™, 380 VAC, 50Hz, 3 Phase

Technical Data

Operation: Manual Chop Cutting

Motor Power: 4.0Hp (3.0kW)

Blade Size and Speed: 10" (254mm) Diameter Blades, 3400 rpm (60Hz), 2830 rpm (50Hz)

Cutting Envelope: Sample Sizes up to 3.75" (95mm)

Controls: Push Button Start and Stop/Emergency Stop; Coolant automatically activates when cutter is started

Other Features: Rock Clamp Assembly, Clean-out Hose, Heavy Duty Interior Lights (14W), Large Viewing Window

Shipping Weight: 200 lbs. (91kg)

Dimensions: 24" W x 17" H x 26.5" D

(600mm x 425mm x 665mm)

Accessories

Blade Options**

11-10074 Continuous Rim Diamond Blade

(for geological materials - yields best quality cuts)

10" x 0.045" x 1.25" (254mm x 1.14mm x 32mm)

Optional Accessories

10-2165-260 Recirculating System, 230 VAC, 60Hz, 3 Phase

10-2165-460 Recirculating System, 460 VAC, 60Hz, 3 Phase

10-2165-250 Recirculating System, 220 VAC, 50Hz, 3 Phase

10-2165-400 Recirculating System, 380 VAC, 50Hz, 3 Phase

16-1497 Recirculating System Cart, with Wheels

Vise Accessories

10-3525 Irregular Specimen Speed Vise Kit (Large)

10-3531 Irregular Specimen Speed Vise Kit (Small)


10-3533 Rock Clamp Assembly



Delta® PetroCut™ Cutter showing Rock Clamp Assembly and Clean-out Hose.



Rock Clamp Assembly with specimen.

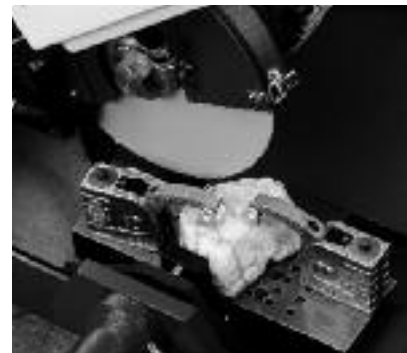

Optional Speed Vise Accessory Kit allows increased flexibility when clamping irregular shapes.

**Blade not included with Delta® PetroCut™

For a complete listing of consumables, please refer to our Consumables Buyer's Guide. Buehler continuously makes product improvements; therefore, technical specifications are subject to change without notice.

© 2010 BUEHLER®, A Division of Illinois Tool Works, Inc. Printed in U.S.A. 25M0106 FN01289 *For metallurgical consumables produced by BUEHLER®



BUEHLER

BUEHLER®, A Division of Illinois Tool Works, Inc.

– Worldwide Headquarters

41 Waukegan Road • P.O. Box 1
Lake Bluff, Illinois 60044-1699 USA
Tel: (847) 295-6500 • Fax: (847) 295-7979
Sales: 1-800-BUEHLER (1-800-283-4537)
www.buehler.com
Email: info@buehler.com



BUEHLER GMBH - European and MESA Headquarters

In der Steele 2 • 40599 Düsseldorf
Postfach 16 03 55 • 40566 Düsseldorf
Telefon: (49) 211 974100 • Telefax: (49) 211 9741079
www.buehler-met.de
Email: info@buehler-met.de

BUEHLER FRANCE

Téléphone: 0800 89 73 71
Télécopie: 0800 88 05 27
www.buehler.fr
Email: info@buehler.fr

BUEHLER UNITED KINGDOM

Telephone: 0800 707 6273
Fax: 0800 707 6274
www.buehler.co.uk
Email: sales@buehler.co.uk

BUEHLER CANADA

10 Carlow Court, Unit #2
Whitby, Ontario L1N 9T7
Telephone: (905) 430-4684
Fax: (905) 430-4647
Sales Telephone: 1-800-268-3593
Email: info@buehler.ca

BUEHLER, ASIA-PACIFIC

5/F Vogue Centre
696 Castle Peak Road
Lai Chi Kok, Kowloon
Hong Kong, SAR, China
Telephone: (852) 2307 0909
Fax: (852) 2307 0233

BUEHLER®

Delta® PetroCut™

GEOLOGICAL CUTTER

- **Affordable, Rugged Manual Cutter Designed for Rocks, Minerals, Nonmetallic and Geological Specimens**
- **Features Universal Clamping with a Rock Clamp Assembly**
- **Produces Serial Cuts Quickly and Consistently**
- **X-Y Table Bed**
- **Specimen Sizes Up To 3.75" (95mm) in Diameter**
- **Uses 10" (254mm) Diamond Blade**
- **Compact Tabletop Design**
- **Heavy-Duty Aluminum Base**
- **Large Viewing Window and Lighted Cutting Chamber**
- **Clean-Out Hose**

