



Delta Medium & Large Service Plans



Protect Your Investment

Labs run more efficiently and consistently with well-maintained equipment.



Reduce Unplanned Downtime

Regular equipment maintenance will reduce the hassle of unexpected downtime and associated costs.



Ensure Consistent Results

Confidently analyze lab results with the equipment maintained by the trusted experts at Buehler.



Preventative Maintenance

Protect your investment and ensure consistent performance, minimize downtime and reduce the costly repairs with Buehler Preventative Maintenance and Calibration Services.

Preventative Maintenance Service Plan Ensures Consistent and Reliable Results

- Annual visit includes Wear Parts Kit and 21 Point Maintenance Check list
- Trained technicians able to perform Preventative Maintenance on multiple machines in one visit
- Unlimited Phone Support





Maintenance Checklist

Our comprehensive 21 point maintenance checklist is thorough, ensuring the Delta Medium/Large is running efficiently.

21 Point Maintenance Checklist

- Check drive & drive pulleys
- Check universal joint
- Check stationary flange
- Check and adjust belt tension
- Lubricate head pivoting pillow blocks
- Check and inspect coolant lines
- Check pivoting pillow block shaft assembly
- Inspect motors and arm assembly for proper operation
- Check for any oxidation or corrosion of the internal motors and internal moving arm assembly
- Clean cutting bay
- Check actuator for wear
- Check and inspect hood & window integrity
- Verify motor and pump rotation for proper operation
- Check over loads for proper amp setting
- Check electronic brake operation and setting
- Check swarf removal system is functioning properly
- Check fume extractor operation
- Check coolant nozzle integrity
- Perform and inspect cuts to ensure machine is working properly
- Check and inspect hood & window integrity
- Check arm bearing housings




Wear Parts Kit

The parts listed below are included with every Preventative Maintenance visit

Delta Medium Wear Parts Kit 102218SRV
 Delta Large Wear Parts Kit 102218SRV

Includes

Chamber Light	Hood Gas Spring
Spindal Lock	Tensioner Bearing
Coolant	Orbital Belt Drive

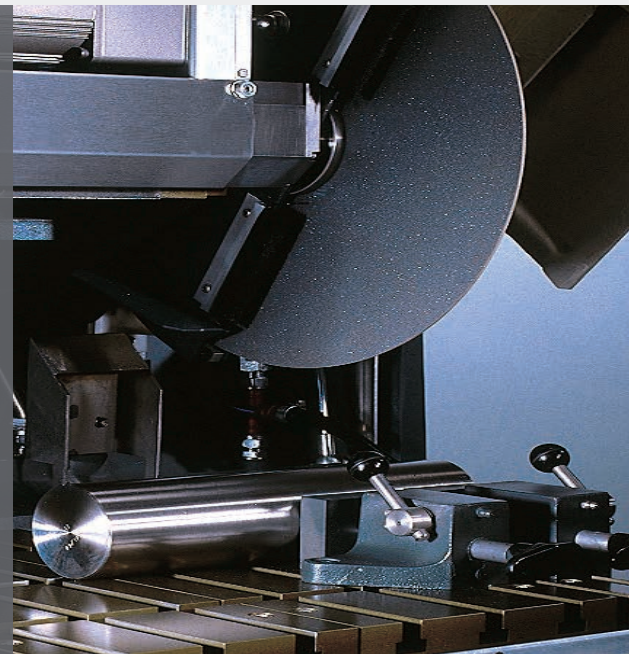

BUEHLER Worldwide Headquarters
North America-South America Offices
 41 Waukegan Road
 Lake Bluff, Illinois 60044-1699 USA
 P: 800 BUEHLER (800-283-4537)
 P: (847) 295-6500
 www.buehler.com | info@buehler.com

European Headquarters
BUEHLER Germany
 info.eu@buehler.com

BUEHLER France
 info.eu@buehler.com

BUEHLER United Kingdom
 info.eu@buehler.com

BUEHLER China
 info.cn@buehler.com

BUEHLER Japan
 info.japan@buehler.com

BUEHLER Asia-Pacific
 info.asia@buehler.com