



EcoMet[®] 30

Manual Single & Twin Service Plans

Protect Your Investment


Labs run more efficiently and consistently with well-maintained equipment.

Reduce Unplanned Downtime

Regular equipment maintenance will reduce the hassle of unexpected downtime and associated costs.

Ensure Consistent Results

Confidently analyze lab results with the equipment maintained by the trusted experts at Buehler.

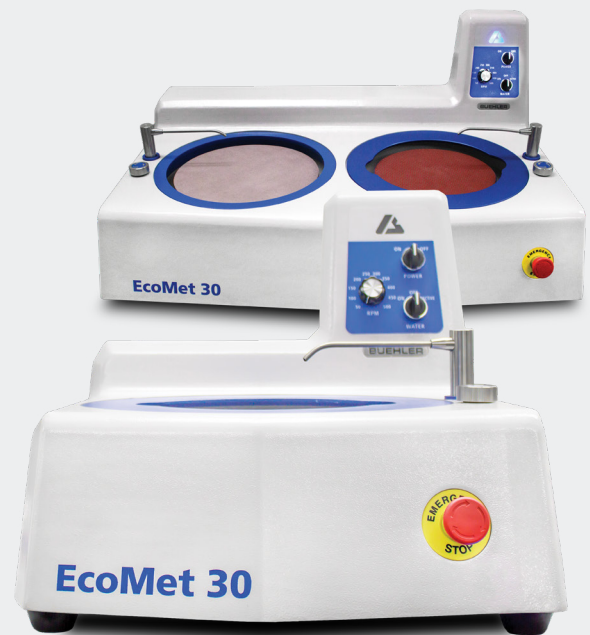


Preventative Maintenance

Protect your investment and ensure consistent performance, minimize downtime and reduce the costly repairs with Buehler Preventative Maintenance and Calibration Services.

Preventative Maintenance Service Plan Ensures Consistent and Reliable Results

- Annual visit includes Wear Parts Kit and 9 Point Maintenance Check list
- Trained technicians able to perform Preventative Maintenance on multiple machines in one visit
- Unlimited Phone Support
- 10% Replacement Part Discount

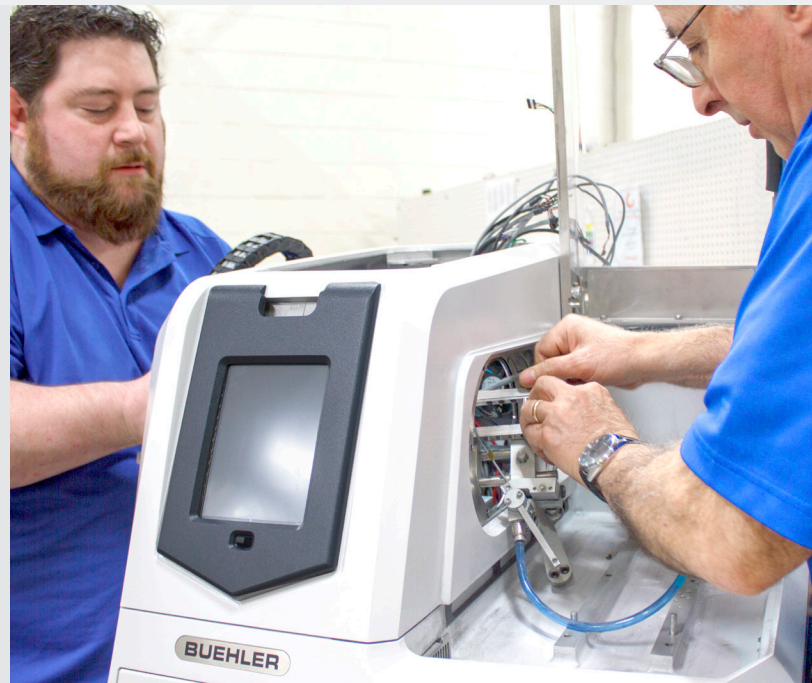



Maintenance Checklist

Our comprehensive 9 point maintenance checklist is thorough, ensuring the EcoMet 30 is running efficiently.

9 Point Maintenance Checklist

- Check that unit components power up properly
- Verify E-stop works properly
- Check display functionality
- Check drain tube
- Check base electrical connections
- Check base motor drive belt
- Check base motor for irregularities and unusual noise
- Check base motor and drive pulley are in proper alignment
- Check water/drain flow and look for any source of leaking



Wear Parts Kit

The parts listed below are among those included with every Preventative Maintenance visit

EcoMet 30 Manual single (4910070) and Twin (4910072) wear parts

Includes

Drain Pipe Assembly
Drive Belt

O-Rings and Oil Seals
Drive Pulley (Twin only)



BUEHLER Worldwide Headquarters
North America-South America Offices
41 Waukegan Road
Lake Bluff, Illinois 60044-1699 USA
P: 800 BUEHLER (800-283-4537)
P: (847) 295-6500
www.buehler.com | info@buehler.com

European Headquarters
BUEHLER Germany
info.eu@buehler.com

BUEHLER France
info.eu@buehler.com

BUEHLER United Kingdom
info.eu@buehler.com

BUEHLER China
info.cn@buehler.com

BUEHLER Japan
info.japan@buehler.com

BUEHLER Asia-Pacific
info.asia@buehler.com