



PlanarMet™ 300 Service Plans



Protect Your Investment

Labs run more efficiently and consistently with well-maintained equipment.



Reduce Unplanned Downtime

Regular equipment maintenance will reduce the hassle of unexpected downtime and associated costs.



Ensure Consistent Results

Confidently analyze lab results with the equipment maintained by the trusted experts at Buehler.




Preventative Maintenance

Protect your investment and ensure consistent performance, minimize downtime and reduce the costly repairs with Buehler Preventative Maintenance and Calibration Services.

Preventative Maintenance Service Plan Ensures Consistent and Reliable Results

- Annual visit includes Wear Parts Kit and 12 Point Maintenance Check list
- Trained technicians able to perform Preventative Maintenance on multiple machines in one visit
- Unlimited Phone Support
- 10% Replacement Part Discount

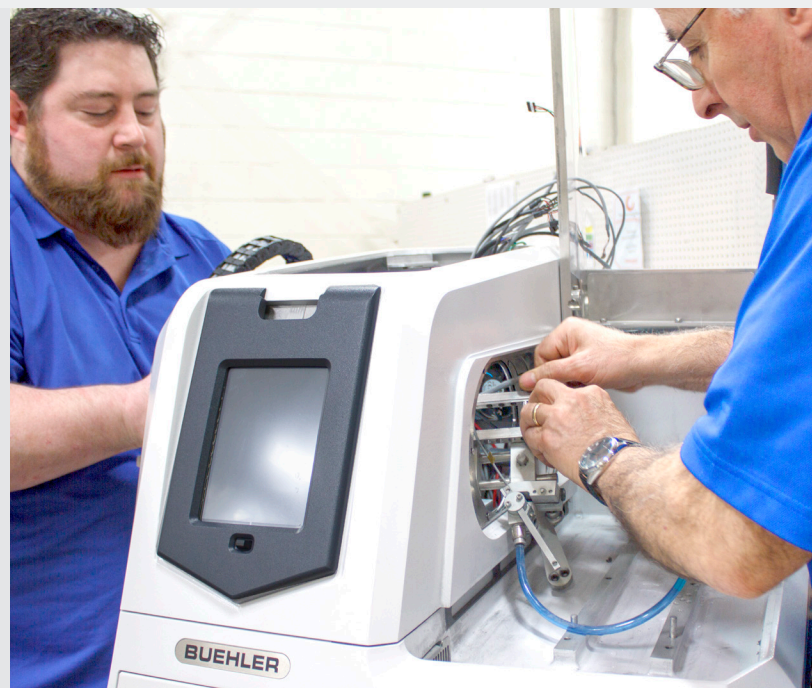



Maintenance Checklist

Our comprehensive 12 Point Maintenance Checklist is thorough, ensuring the PlanarMet 300 is running efficiently.

12 Point Maintenance Checklist

- | | |
|---------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------|
| <input checked="" type="checkbox"/> Check tension and condition of platen, power head and dressing belt | <input checked="" type="checkbox"/> Properly lubricate power head components and dressing system |
| <input checked="" type="checkbox"/> Inspect bearing carrier seals | <input checked="" type="checkbox"/> Check to see if grease is leaking from bearing carrier |
| <input checked="" type="checkbox"/> Observe machine for proper operation | <input checked="" type="checkbox"/> Inspect splash cover manifold for debris build-up |
| <input checked="" type="checkbox"/> Inspect splash cover gasket | <input checked="" type="checkbox"/> Inspect dressing boot for tears or leaks |
| <input checked="" type="checkbox"/> Clean bowl and drain tube for proper drainage | <input checked="" type="checkbox"/> Adjust power head safety zone for proper operation |
| <input checked="" type="checkbox"/> Clean chuck sleeve (moves with ease) | <input checked="" type="checkbox"/> Adjust platen motor brake, if slow to stop |




Wear Parts Kit

The parts listed below are included with every Preventative Maintenance visit

PlanarMet 300 Wear Parts Kit 491000SRV

Includes

Splash Cover Gasket	Idler Pulley
Platen Carrier V-Ring Seal	Dressing Stick
Platen Carrier Inner Ring Seal	Coolant
Drive Belt	Coolant O-Ring

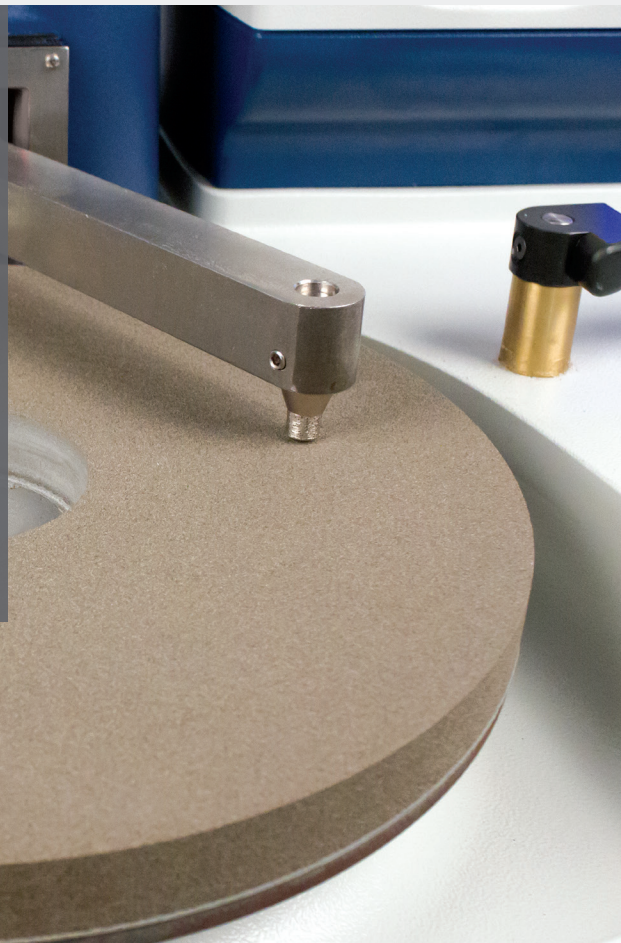

BUEHLER Worldwide Headquarters
North America-South America Offices
41 Waukegan Road
Lake Bluff, Illinois 60044-1699 USA
P: 800 BUEHLER (800-283-4537)
P: (847) 295-6500
www.buehler.com | info@buehler.com

European Headquarters
BUEHLER Germany
info.eu@buehler.com

BUEHLER France
info.eu@buehler.com

BUEHLER United Kingdom
info.eu@buehler.com

BUEHLER China
info.cn@buehler.com

BUEHLER Japan
info.japan@buehler.com

BUEHLER Asia-Pacific
info.asia@buehler.com